
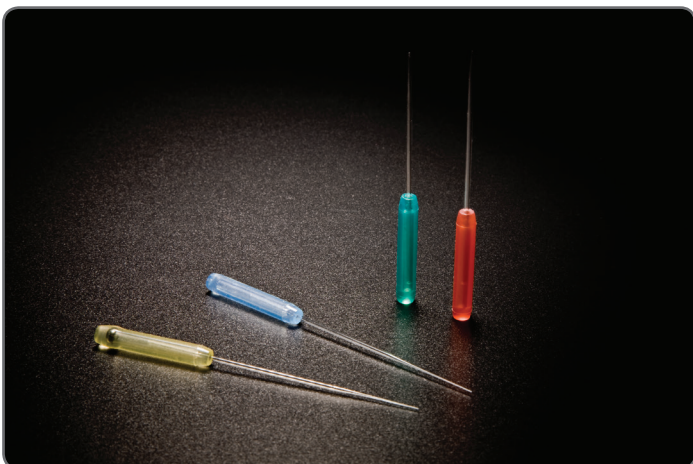
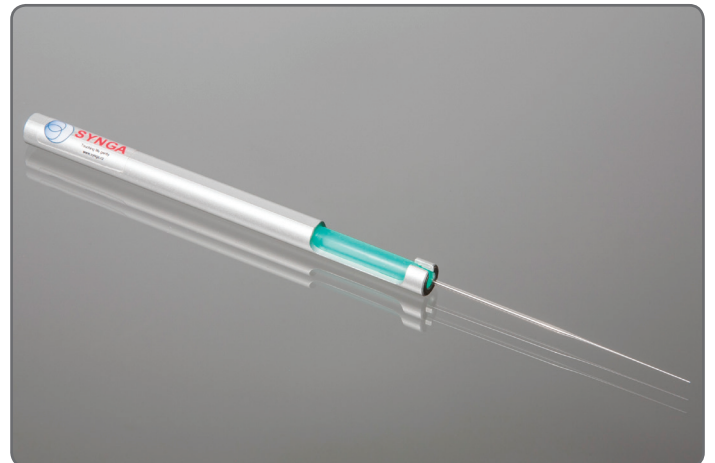
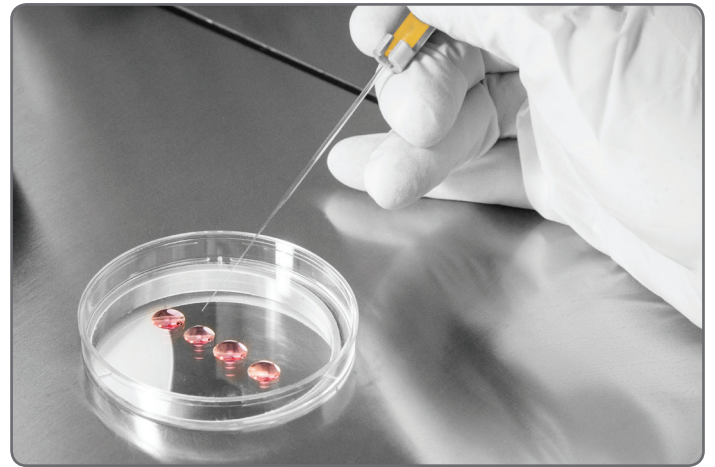


# Synga Premounted

---

**Denudation and Manipulation System  
for Streamlined Routines in IVF Labs**

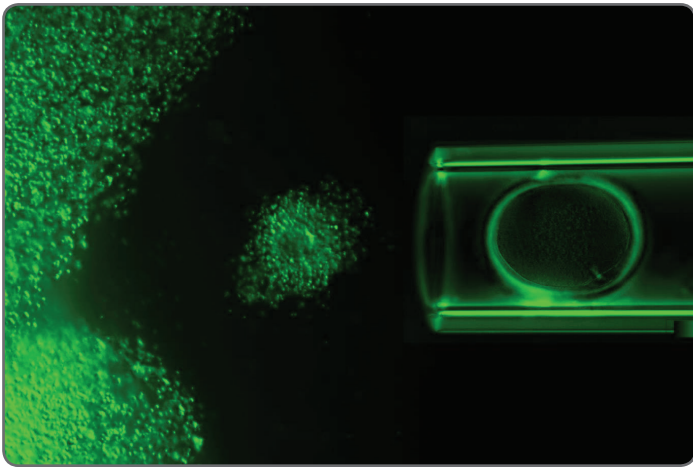

## Synga Premounted

Synga Premounted is a pipette and bulb system for oocyte denudation or gamete and embryo handling, designed for use with a handle accessory (penholder).

Each pipette is fitted with an aspiration bulb, creating a reliable seal, which is not always guaranteed with customer assembled systems. This time-saving combined unit enables a faster set-up. Colour coding of bulbs allows straightforward differentiation of diameters. High quality borosilicate glass minimizes the risk of oil attaching to the inner surface – no oil clogging.

Synga Premounted pipettes are supplied sterile in individually sealed blister packs. Each pipette remains sterile until its blister card cover is peeled off open. The autoclavable Synga penholder (handle), supplied separately, is designed to offer precision, ease of handling and a comfortable working position.

To achieve maximum quality, each pipette is inspected several times during manufacture. Visual check for absence of glass fragments is followed by automated measurement of diameter to keep it within specified tolerance.

### Two types of bulbs - Soft or VI

The pipettes come with specially designed sterile silicone bulbs. The bulb has a smaller volume than the pipette to eliminate the risk of oocyte trapping. VI bulbs are thinner and softer to allow for individual needs.


### Two lengths of pipettes - S or L

Denudation diameters of pipettes are produced either with functional part 14 mm long, marked S in reference number, or with functional part 29 mm long, marked L.

## Synga Accessories



Penholders

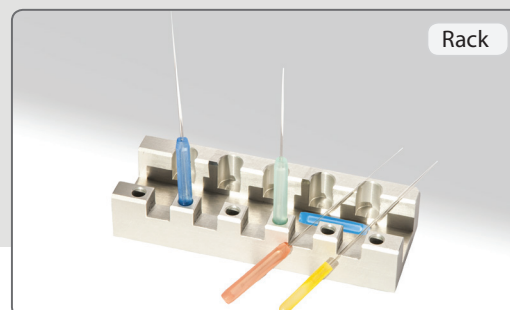


Rack

Description	Inner diameter	Tip	Colour
Denudation SG Pipette	Ø 130 µm	sharp	red ●
Denudation SG Pipette	Ø 140 µm	sharp	green ●
Denudation SG Pipette	Ø 150 µm	sharp	white ○
Vitrification SG Pipette	Ø 150 µm	polished	white ○
Manipulation SG Pipette	Ø 180 µm	polished	yellow ●
Blastocyst SG Pipette	Ø 250 µm	polished	blue ●

**Products and innovations from experienced embryologists**

We proudly serve life-focused laboratories.

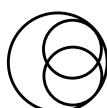

**SYNGA**  
TOUCHING LIFE GENTLY